

Shoulder Cryo/Cuff™ IC Patient Information

Patient: _____

Doctor: _____

What is the Cryo/Cuff IC?

The Aircast Cryo/Cuff IC (Intermittent Compression) is an integrated approach to cold therapy that combines cold and focal compression to better manage pain and reduce swelling. Intermittent cold and focal compression are delivered using a motorized delivery system.

How is the Cryo/Cuff IC used?


STEP 1 - Before applying the Cryo/Cuff, prepare the cooler (Figure 1)

1. Connect the blue tube to cooler.
2. Add cold water to line inside cooler.
3. Fill with ice.
4. Lay insulation disk on top of ice. Attach cooler lid snugly.
5. Allow five minutes with occasional shaking to chill water.



STEP 2 - Always apply EMPTY Shoulder Cryo/Cuff (Figure 2) (The first time, have another person help)

- Use only Aircast CryoCuffs and ArcticFlow Cuffs with Aircast Cryo/Cuff IC Unit.
- Snugly secure strap around torso in front and back. Adjust underarm strap (first time only). Secure front and back elastic straps.
- Carefully and completely read instructions supplied with the specific Cryo/Cuff prior to application. Please follow the instructions below to continue cooler preparation.

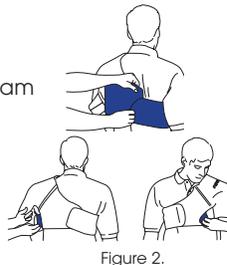

STEP 3 - Fill and pressurize Cryo/Cuff (Figure 3)

- To fill cuff-connect blue tube to Cryo/Cuff. Raise cooler no more than 15" (38 cm) above the cuff for approximately 30 seconds while Cryo/Cuff fills.
- To correctly position cooler-place cooler on stable surface. The cooler base should be even with the Cryo/Cuff - NOTE: To AVOID excessive pressure during use, DO NOT raise the cooler higher than 15" (38 cm) above the Cryo/Cuff.
- To turn cooler on-plug power supply into wall outlet and connect cord to lid.



⚠ CAUTION

- DO NOT use an elastic wrap with the Cryo/Cuff IC.
- Dressings used under the Cryo/Cuff should be applied lightly.
- When filling the Cryo/Cuff IC, DO NOT raise the cooler higher than 15" (38 cm) above the Cryo/Cuff avoid excessive pressure.
- Reduce pressure with any sense of discomfort, numbness or tingling of the limb.
- EMPTY the Cryo/Cuff after each use, and prior to refilling.
- Patient skin condition checks must occur every 1-2 hours on a routine basis.

To re-fill cooler:

When ice has melted, unplug the cooler.

IMPORTANT: EMPTY the Cryo/Cuff COMPLETELY while still attached to the cooler.

When Cryo/Cuff is empty then empty the cooler. Repeat steps 1-5 in **STEP 1** to prepare the cooler for refilling.

To clean cooler:

Use mild detergent to clean cooler and lid. DO NOT submerge lid under water.

To troubleshoot:

1. Make sure unit has ice and water filled to the line indicated inside the cooler.
2. Check power connections at wall and at the unit.
3. Make sure to fill Cryo/Cuff prior to operating unit as instructions indicate.
4. Check Cryo/Cuff connection; make sure to listen for a click when connecting the hose to the Cryo/Cuff.
5. Make sure unit and Cryo/Cuff are level for optimal performance.
6. This device may be susceptible to magnetic or electromagnetic interference. If this occurs, relocate or reposition the device.
7. Increase the separation between the equipment.
8. Plug the power supply into an outlet on a circuit different than other devices.
9. If unit motor is still not functioning-unit can be used as a gravity fed device (see Figure 3) for positioning.



DJOglobal.com

FOR PRODUCT SUPPORT CALL
1-888-405-3251

AIRCAST®

DJO, LLC.
1430 Decision Street, Vista, CA 92081-8553 U.S.A.
Phone 760.727.1280 Toll Free 800.526.8785 Fax 800.457.4221
aircast.com