
# X-ACT ROM ELBOW

- Codo X-Act ROM
- X-Act ROM Ellenbogen
- Gomito X-Act ROM
- Attelle de coude X-Act ROM Elbow
- X-act ROM elleboog
- X-Act ROM Elbow
- X-Act ROM Elbow
- X-Act ROM Elbow
- X-Act ROM Elbow
- X-Act ROM -kyynärpäätuki
- X-Act ROM Elbow
- Cotovelo X-Act ROM
- X-act ROM albue
- X-Act ROM Elbow
- X-Act ROM エルボー

ENGLISH .....	2
ESPAÑOL .....	4
DEUTSCH .....	6
ITALIANO .....	8
FRANÇAIS .....	10
NEDERLANDS .....	12
SLOVENČINA .....	14
DANSK .....	16
SVENSKA .....	18
SUOMI .....	20
ČEŠTINA .....	22
PORTUGUÊS .....	24
NORSK .....	26
简体中文 .....	28
日本語 .....	30

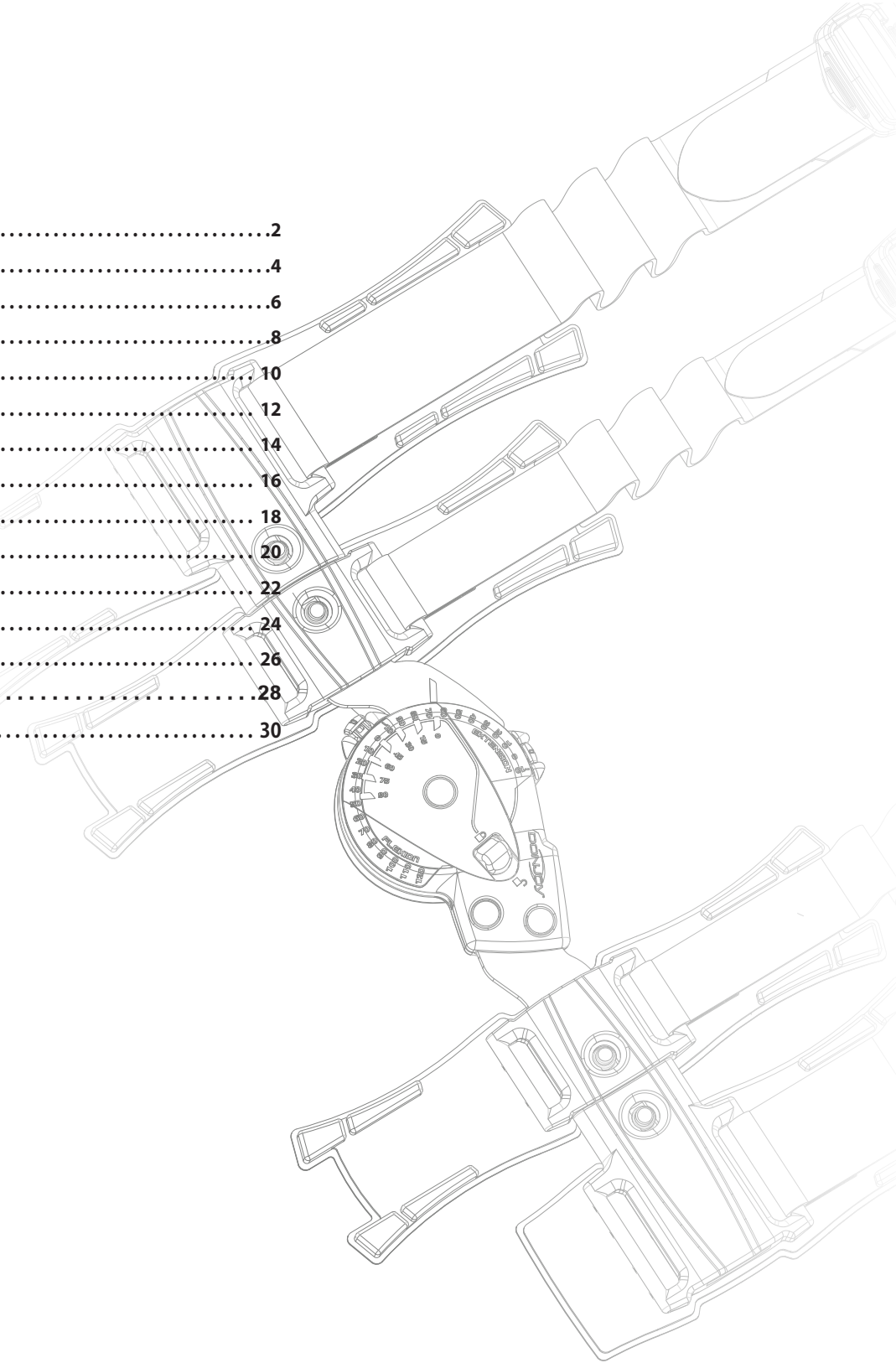DJO, LLC  
A DJO Global Company  
1430 Decision Street  
Vista, CA 92081-8553 • USA

T 760.727.1280 • F 760.734.5608  
DJOglobal.com

13-6811\_A

**BEFORE USING THE DEVICE, PLEASE READ THE FOLLOWING INSTRUCTIONS COMPLETELY AND CAREFULLY. CORRECT APPLICATION IS VITAL TO THE PROPER FUNCTIONING OF THE DEVICE.**

**INTENDED USE/INDICATIONS:** Provides fixed or range of motion control post-injury or post-operative procedure to the elbow; stable fractures; biceps repair; and mild to severe soft tissue injuries.

**CONTRAINDICATIONS:** None

**WARNINGS AND PRECAUTIONS:** If you experience any pain, swelling, sensation changes, or any unusual reactions while using this product, consult your medical professional immediately.

**APPLICATION INFORMATION:**

1. **Unlock Hinge.** The X-Act Rom Elbow is packaged in a locked position. To unlock brace, slide Red Quick-Lock button away from hinge center.
2. **Position Brace.** Place the mid-point of hinge over lateral epicondyle (outside of elbow) **Note: Red Quick Lock Button should be below Elbow when brace is applied.**
3. **Telefit Arm Lengths.** While hinged is aligned with epicondyle (mid-point of elbow) adjust position of all four cuffs. Depress the slider button and position both forearm cuffs then both bicep cuffs to desired location. **Note:** Make sure proximal Bicep cuff extends to most comfortable position, but does not impinge Axilla (arm pit) and distal Forearm cuff is positioned above styloid process (bump on lateral side of wrist). **DO NOT attempt to remove Sliders from hinge bar.**
- 4a. **Apply Brace.** Unfasten the four Quick-Clip buckles on sliders by rotating the release end upward until clip is ejected off of slider. Place arm into brace.
- 4b. Starting with the cuffs closest to elbow, form cuffs around forearm and bicep.
5. **Secure Straps.** Starting at the wrist, click each of the Quick-Clip buckles onto the slider bar, making sure an “audible click” is achieved. Cinch Straps for snug fit.  
*Trim Straps as necessary.* Remove alligator strap tab; trim strap; replace strap tab. **NOTE:** Bicep Strap pads are packed with brace for comfort and anti-migration.
6. **Range of Motion Hinge Setting.** X-Act Elbow brace allows for full range of motion in 10 degree increments from -10 to 120 degrees of flexion. To adjust settings, depress the Flexion & Extension buttons and move to prescribed range of motion (ROM).
7. **Quick-Lock.** To prevent any range of motion, slide the red Quick-Lock button to the locked position. X-Act Elbow brace may be locked in 0, 15, 30, 45, 60, 75, or 90 degrees of flexion.

8. **Patient Compliance.** If required, lock ties are included with each brace. Loop a tie through hole in both the extension and flexion buttons to maintain prescribed range of motion settings.
9. **Removal & Re-Application of X-Act Rom Elbow Brace.** Slightly loosen straps on each cuff and disengage clip on sliders by rotating clip upward until clip is ejected off of slider. Loosen cuffs. **Re-apply brace easily with one hand. Align hinge at Elbow.** Form cuffs around arm for snug fit. Click buckles onto sliders. Tighten straps for proper fit.

**APPLYING SHOULDER SLING:**

10. Loosen both straps on Shoulder Sling, Place shoulder strap over head so strap sits on both shoulders, like a collar.
11. Place uninjured arm through Shoulder Sling, making sure one strap is across back. Shoulder strap should be over shoulder of injured arm.
12. Attach Shoulder Sling to brace by feeding wrist strap through buckle of Shoulder Sling. Wrap strap around wrist cuff and secure.
13. Move contact closure strap ends on Shoulder Sling to raise arm to comfortable position. Strap may be trimmed as needed. Position shoulder pad on neck area as needed for comfort.

**To remove sling** - Use Quick Release buckles for easy removal and re-application of product

**CLEANING INSTRUCTIONS:**

Remove and hand wash preformed cuff liners and brace straps in cold water with mild detergent. Rinse thoroughly. AIR DRY only, do not heat dry. Regular cleaning of the brace is recommended.

**WARRANTY:** DJO, LLC will repair or replace all or part of the unit and its accessories for material or workmanship defects for a period of six months from the date of sale.

**FOR SINGLE PATIENT USE.  
RX ONLY**

NOTICE: WHILE EVERY EFFORT HAS BEEN MADE IN STATE-OF-THE-ART TECHNIQUES TO OBTAIN THE MAXIMUM COMPATIBILITY OF FUNCTION, STRENGTH, DURABILITY AND COMFORT, THIS DEVICE IS ONLY ONE ELEMENT IN THE OVERALL TREATMENT PROGRAM ADMINISTERED BY A MEDICAL PROFESSIONAL. THERE IS NO GUARANTEE THAT INJURY WILL BE PREVENTED THROUGH THE USE OF THIS PRODUCT.

CAUTION: FEDERAL LAW (U.S.A.) RESTRICTS THIS DEVICE TO SALE BY OR ON THE ORDER OF A LICENSED HEALTH CARE PROFESSIONAL.

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