enovis

Active, adaptive compression for subtalar fusion in an easy-to-insert design.

DYNANAIL HYBRIDTM

FUSION SYSTEM

FEATURES & BENEFITS:

- Hybrid screw/nail design combines dynamic, compressive power with simple, intuitive insertion.
- Internal NiTiNOL element offers between 2 5 mm of active compression in response to bone settling or resorption, PLUS manual intraoperative compression.
- 16 mm proximal talar threads allow for easy insertion, similar to a screw.
- Distal transverse screw provides additional stability and prevents device migration.
- NiTiNOL Compressive Element comes pre-stretched and pre-loaded on a disposable Nail Guide for enhanced efficiency.



DYNANGE.


UNIVERSAL INSTRUMENTATION

Shares the same instrumentation as the DynaNail Mini® to provide ease-of-use and intra-operative flexibility.

THE ULTIMATE IN MAINTAINED **ACTIVE COMPRESSION**

Unlike traditional screws that lose compression post-operatively, the DynaNail Hybrid™ maintains compression up to 5 mm of bone resorption or settling depending on implant length.1

COMPRESSION PERFORMANCE SUMMARY

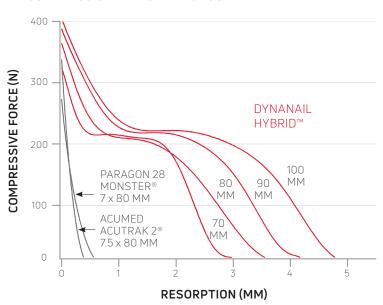

DYNANAIL HYBRID™ IMPLANTS

DESCRIPTION	PART #
DYNANAIL HYBRID, 7 MM X 60 MM	2601-00-7060
DYNANAIL HYBRID, 7 MM X 70 MM	2601-00-7070
DYNANAIL HYBRID, 7 MM X 80 MM	2601-00-7080
DYNANAIL HYBRID, 7 MM X 90 MM	2601-00-7090
DYNANAIL HYBRID, 7 MM X 100 MM	2601-00-7100
HEADLESS SCREW, 4 MM X 22 MM	2600-04-4322
HEADLESS SCREW, 4 MM X 24 MM	2600-04-4324
HEADLESS SCREW, 4 MM X 26 MM	2600-04-4326
HEADLESS SCREW, 4 MM X 28 MM	2600-04-4328
HEADLESS SCREW, 4 MM X 30 MM	2600-04-4330
HEADLESS SCREW, 4 MM X 32 MM	2600-04-4332
HEADLESS SCREW, 4 MM X 34 MM	2600-04-4334
HEADLESS SCREW, 4 MM X 36 MM	2600-04-4336
HEADLESS SCREW, 4 MM X 38 MM	2600-04-4338
HEADLESS SCREW, 4 MM X 40 MM	2600-04-4340

DYNANAIL HYBRID™ SINGLE-USE INSTRUMENTS

DESCRIPTION	PART #
STEINMANN PIN, 2 MM X 9"	2200-19-0020
2.5 MM X 6" DRILL	2201-09-0025
TRANSVERSE SCREW DRILL, 4 MM, STEPPED	2900-03-0400
GUIDEWIRE, 2.4 MM X 229 MM	2900-04-0229
CANNULATED STEP DRILL, 7.5 MM	2901-01-0075

CHECK OUT MORE OF THE DYNANAIL FAMILY

WITH UNIVERSAL INSTRUMENTATION



DYNANAIL HYBRID™ **FUSION SYSTEM**

FUSION SYSTEM

T 800.495.2919 F 877.778.3864

1575 Northside Drive NW | Suite 440 | Atlanta, GA 30318 | U.S.A. enovis.com/foot-and-ankle

1. Data on File. Laboratory testing does not necessarily indicate clinical performance.

Individual results may vary. Neither Medshape, Inc. nor any of the Enovis companies dispense medical advice. The contents of this document do not constitute medical, legal, or any other type of professional advice. Rather, please consult your healthcare professional for information on the courses of treatment, if any, which may be appropriate for you.