

ARSENAL FOOT PLATING SYSTEM™

PLATING SYSTEM

The Arsenal Foot Plating System™ equips surgeons with cutting-edge technology that redefines the benchmark of plating. Created from Type II Anodized Titanium Alloy, this comprehensive plating system is designed to effectively manage fusions, fractures, and osteotomies across the forefoot and midfoot.

DESIGN FEATURES

- 60° locking cone of angulation surpasses competitors by offering up to 30° in any direction.
- Lock and relock locking screws up to five times, a design feature aimed toward minimizing implant waste.
- Dual lead 2.7mm and 3.5mm locking and non-locking screws are designed to save valuable OR time.
- Full offering for MTP, Lapidus, Lisfranc, and Utility plates meet variety of patient anatomies and indications.
- Made from Type II Anodization, Titanium Alloy for optimized thickness and plate contouring.
- Scalloped plate design facilitates easy contouring to the patient's anatomy while minimizing plate footprint.

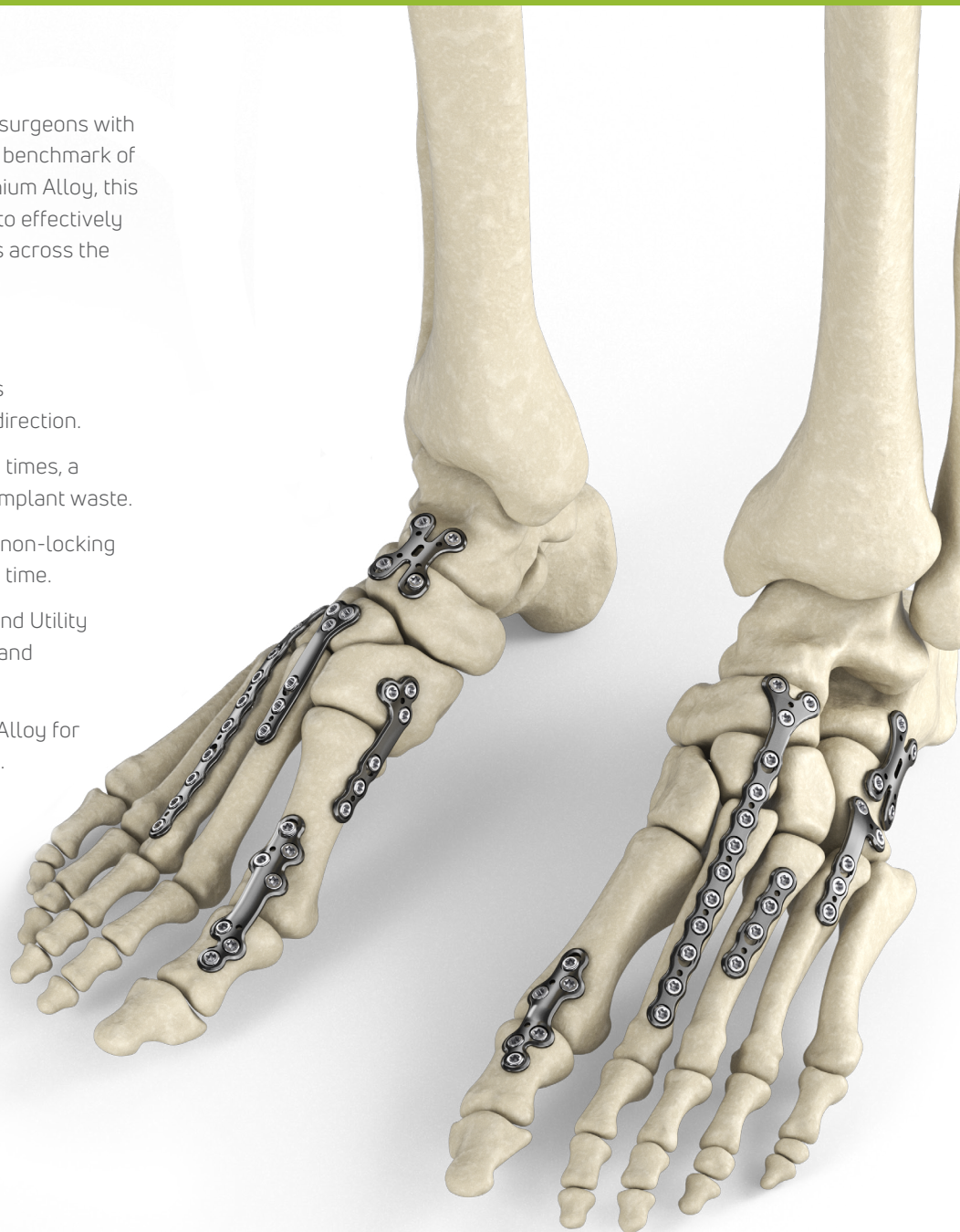


PLATE FAMILIES



MTP PLATES – 16 OFFERINGS

- 0° and 5° options
- Petite, standard, long, and revision plates



LAPIDUS PLATES – 6 OFFERINGS

- Petite, short, and long plates



LISFRANC PLATES – 6 OFFERINGS

- 0°, 10°, and 20° options
- Standard and long plates



UTILITY PLATES – 19 OFFERINGS

- Straight, L, X, and Y Plates

ARSENAL FOOT PLATING SYSTEM™ ORDERING INFORMATION

PART #	DESCRIPTION
300-82-011	PETITE MTP PLATE, L, 5°
300-82-012	PETITE MTP PLATE, R, 5°
300-82-013	PETITE MTP PLATE, L, 0°
300-82-014	PETITE MTP PLATE, R, 0°
300-82-015	STANDARD MTP PLATE, L, 5°
300-82-016	STANDARD MTP PLATE, R, 5°
300-82-017	STANDARD MTP PLATE, L, 0°
300-82-018	STANDARD MTP PLATE, R, 0°
300-82-019	LONG MTP PLATE, L, 5°
300-82-020	LONG MTP PLATE, R, 5°
300-82-021	LONG MTP PLATE, L, 0°
300-82-022	LONG MTP PLATE, R, 0°
300-82-023	REVISION MTP PLATE, L, 5°
300-82-024	REVISION MTP PLATE, R, 5°
300-82-025	REVISION MTP PLATE, L, 0°
300-82-026	REVISION MTP PLATE, R, 0°
300-83-007	PETITE LAPIDUS PLATE, LEFT
300-83-008	PETITE LAPIDUS PLATE, RIGHT
300-83-009	SHORT LAPIDUS PLATE, LEFT
300-83-010	SHORT LAPIDUS PLATE, RIGHT
300-83-011	LONG LAPIDUS PLATE, LEFT
300-83-012	LONG LAPIDUS PLATE, RIGHT
300-80-007	2-HOLE STRAIGHT PLATE
300-80-008	4-HOLE STRAIGHT PLATE

PART #	DESCRIPTION
300-80-009	5-HOLE STRAIGHT PLATE
300-80-010	6-HOLE STRAIGHT PLATE
300-80-011	8-HOLE STRAIGHT PLATE
300-80-012	10-HOLE STRAIGHT PLATE
300-89-001	4-HOLE Y PLATE
300-89-002	5-HOLE Y-PLATE
300-89-003	7-HOLE Y-PLATE
300-89-004	9-HOLE Y-PLATE
300-89-005	12-HOLE Y-PLATE
300-88-003	SHORT 4-HOLE L-PLATE, LEFT
300-88-004	SHORT 4-HOLE L-PLATE, RIGHT
300-88-005	4-HOLE L-PLATE, LEFT
300-88-006	4-HOLE L-PLATE, RIGHT
300-88-007	5-HOLE L-PLATE, LEFT
300-88-008	5-HOLE L-PLATE, RIGHT
300-88-009	STANDARD X-PLATE
300-88-010	LARGE X-PLATE
300-85-009	STANDARD LISFRANC, 20°, PLATE
300-85-010	STANDARD LISFRANC, 10°, PLATE
300-85-011	STANDARD LISFRANC, 0°, PLATE
300-85-012	LONG LISFRANC, 20°, PLATE
300-85-013	LONG LISFRANC, 10°, PLATE
300-85-014	LONG LISFRANC, 0°, PLATE

enovis

T 800.495.2919 F 877.778.3864

Trilliant, LLC
727 North Shepherd Dr | Suite 100 | Houston, TX 77007 | U.S.A.
djglobal.com/foot-and-ankle

Copyright © 2024 Enovis Foot and Ankle
MK-10533 Rev 0

Individual results may vary. Neither Trilliant, LLC nor any of the Enovis companies dispense medical advice. The contents of this document do not constitute medical, legal, or any other type of professional advice. Rather, please consult your healthcare professional for information on the courses of treatment, if any, which may be appropriate for you.