



#1 in Knee Bracing

CCMI Mark III



DJO Global

Vista Manufacturing Center

3151 Scott Street | Vista, CA 92081-8553 USA

Tel: (800) 336-6569 | www.djoglobal.com

Fax To: (877) 595-8623

Email To: A22.help@djoglobal.com

BILL TO:

Account Name: \_\_\_\_\_

Account #: \_\_\_\_\_ P.O. #: \_\_\_\_\_

Contact: \_\_\_\_\_ Phone: \_\_\_\_\_

Address: \_\_\_\_\_

City: \_\_\_\_\_ State: \_\_\_\_\_ Zip: \_\_\_\_\_

SHIPPING INFORMATION:  Same as "BILL TO"

Name: \_\_\_\_\_

Address: \_\_\_\_\_

City: \_\_\_\_\_ State: \_\_\_\_\_ Zip: \_\_\_\_\_

Red (overnight)  Blue (2 day air)

3 Day  Ground  Other \_\_\_\_\_

PRODUCT SELECTION: CCMI Mark III A22

FourcePoint™ (ACL only)

FEATURES for A22

Installed Extension Stop: 0° 10° 20° 30° 40° (circle one)

Installed Flexion Stop: 45° 60° 75° 90° (circle one)

Neoprene Suspension Strap Kit

ACCESSORIES (extra charge):

A22 Sleeve ( S  M  L)

Full Extension Lockout Stop  Armor Impact Guard

Lycra® Suspension Sleeve  Pinch Guard (Calf)

Neoprene Suspension Sleeve  Atrophy Pad Kit

Sports Brace Cover  Air Condyle Kit

A22 Gel Condyle Pads

FOR DJO USE ONLY:

Order #: \_\_\_\_\_

Brace #: \_\_\_\_\_

PATIENT INFORMATION:

Name: \_\_\_\_\_

Age: \_\_\_\_\_ Height: \_\_\_\_\_ Weight: \_\_\_\_\_

Knee Measurement:  Right  Left  Reverse

Measured By: \_\_\_\_\_


New Brace

WARRANTY:

Check One

Collision Sports (American Football) - 6 months

Standard warranty -1 year

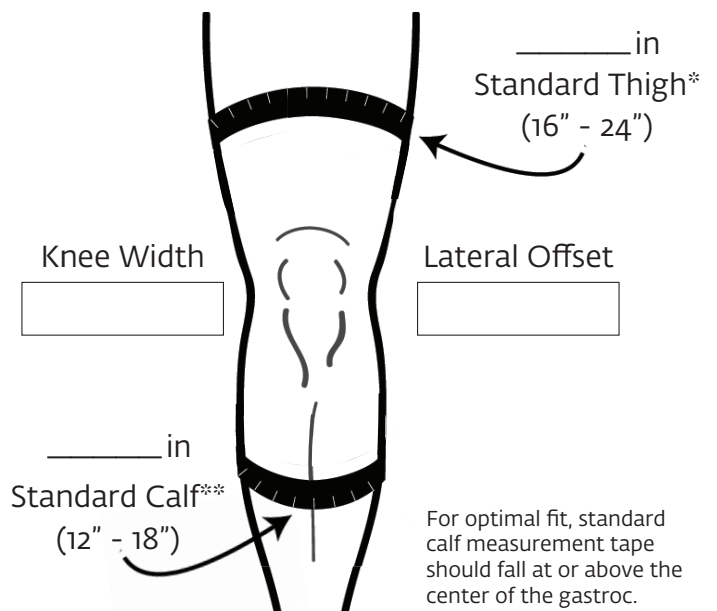

\*Only available in standard thigh

\*\*Record only standard calf measurements (actual calf cuff length is a hybrid between short and standard)


For further instructions, please refer to the CCMI Mark III instruction video and instruction card included with the instrument.

## CAUTION!

FEDERAL LAW (U.S.A.) RESTRICTS THIS DEVICE TO SALE BY THE ORDER OF A LICENSED HEALTH CARE PROFESSIONAL

### Custom Contour Measuring Instrument - Mark III Use Instructions

NOTICE: For optimum long term fit of the DEFIANCE® or DEFIANCE® III and A22, measurements should be taken when there is no edema or atrophy present.

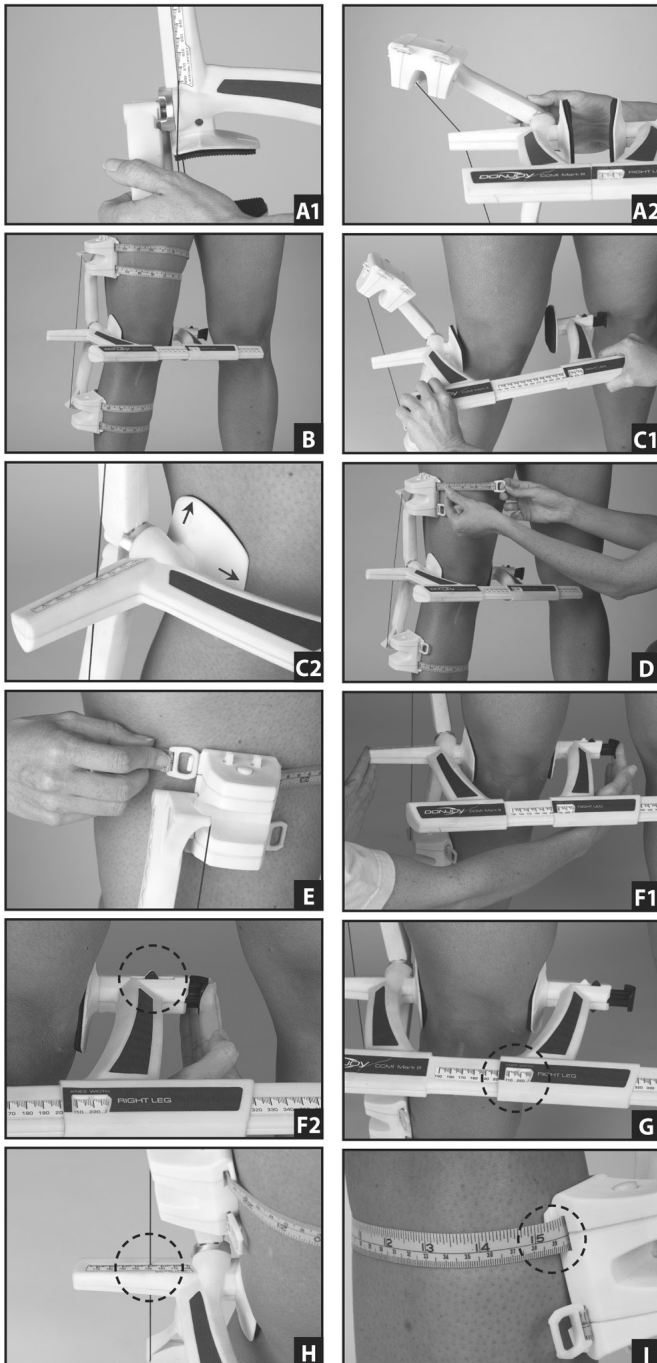

#### Application Information for A22:

1. To unstow, distract and rotate swing arms as shown in figures (A1 and A2).
2. Have the patient stand feet shoulder width apart and knees flexed slightly. Do not allow the knees to flex beyond 15 degrees. Apply the instrument directly to the bare leg. Note: Taking measurements from the uninvolved leg is not recommended because of differences in leg dimensions due to dominance of one leg over the other. See figure B.
3. Distract knee width caliber and position on the patient's knee. Arrows on condyle shells should align with the top of the knee cap and slightly posterior to midline. See figures (C1 and C2).
4. For standard thigh and calf, extend tape measure marked in Inches located distal to knee. Wrap around leg and secure to hook on tape module. Press tape release button to pull tape snug to leg; ensure it is parallel to floor. **“For A22, use only standard thigh and standard calf measurements.”** See figures (D and E).
5. To establish proper knee compression for knee width measurement, support lateral side of instrument and press load plunger on medial condyle until triangles are aligned, as shown in figures (F1 and F2).
6. Record knee width on order form. See figure (G.)
7. Record lateral offset on order form. See figure (H).
8. Record thigh and calf circumference on order form. See figure (I).